


Application Guide

XLR8/XLR8 Plus

GENADYNE

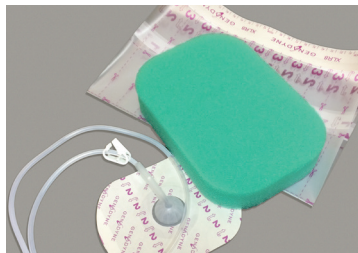

Step 1

Select appropriate size dressing kit for wound. Kit contains foam, port pad, and 2 drapes. (Small kit has 1 drape only.)



Step 2

Clean the wound according to protocol. Protect wound edges with drape.



Step 3

Cut foam to fit accurately in wound bed. Avoid overlapping onto intact skin. Avoid cutting foam over wound. Avoid over-packing, as this will inhibit filtration of drainage.



Step 4

Cover foam with drape by separating bottom layers 1 and 2 to expose the adhesive drape. Position and place drape over foam. Remove stabilization layers 1 and 2, then top with layer 3. Remove white ruler tab.



Step 5

Cut a hole in drape the approximate size of a quarter. Remove paper backing 1 from bottom of the Port Pad. Position Port Pad opening over the hole that was created. Remove backing 2 from the top of Port Pad. Flip handler tabs over on each side and peel down to remove.



Step 6

Select canister size and attach to the back of the pump. Connect the port pad to the canister by securing the purple luer lock.



Step 7

Initiate therapy by pressing the universal on/off button. Adjust pressure as prescribed using the up/down arrows.



Step 8

A secure seal is identified when the foam draws down and the target pressure is reached.



Questions? Call: 732-367-5533 | or email: info@handrhealthcare.com

www.handrhealthcare.com