

THE GENADYNE All in One HYBRID DRESSING

For use with XLR8+ and UNO+ NPWT Pump



A convenient dressing system that is EASY to use, QUICK to apply, COMFORTABLE for the patient, and CLINICALLY effective in treating various types of wounds.

Easy and Quick Simply Peel off the silicone backing and place over the wound and no cutting required

Comfortable The Silicone layer provides comfort and flexibility

Positionable The Silicone layer allows for easy access to the wound with ability to reapply

Cost Savings if used alone may stay on the patient up to 7 days

Insurance Coverage Billable for homecare under code A6550

Wound Types Post-Operative incision management (wounds with 1-2cm depth, wounds over 2cm may pack Genadyne's fully reticulated polyurethane foam below the hybrid dressing). Plastic and Reconstructive, orthopedic, general surgery, podiatry, and vascular


Two Linear dressing options and Two wider rectangular dressing options both for one step application over a multitude of wounds (burns, vascular, diabetic) and skin grafts

10cm x 30cm
10cm x 50cm
15cm x 20cm
15cm x 30cm



UNO HYBRID DRESSING

The all in one dressing system for use with Genadyne's XLR8+ and Uno+

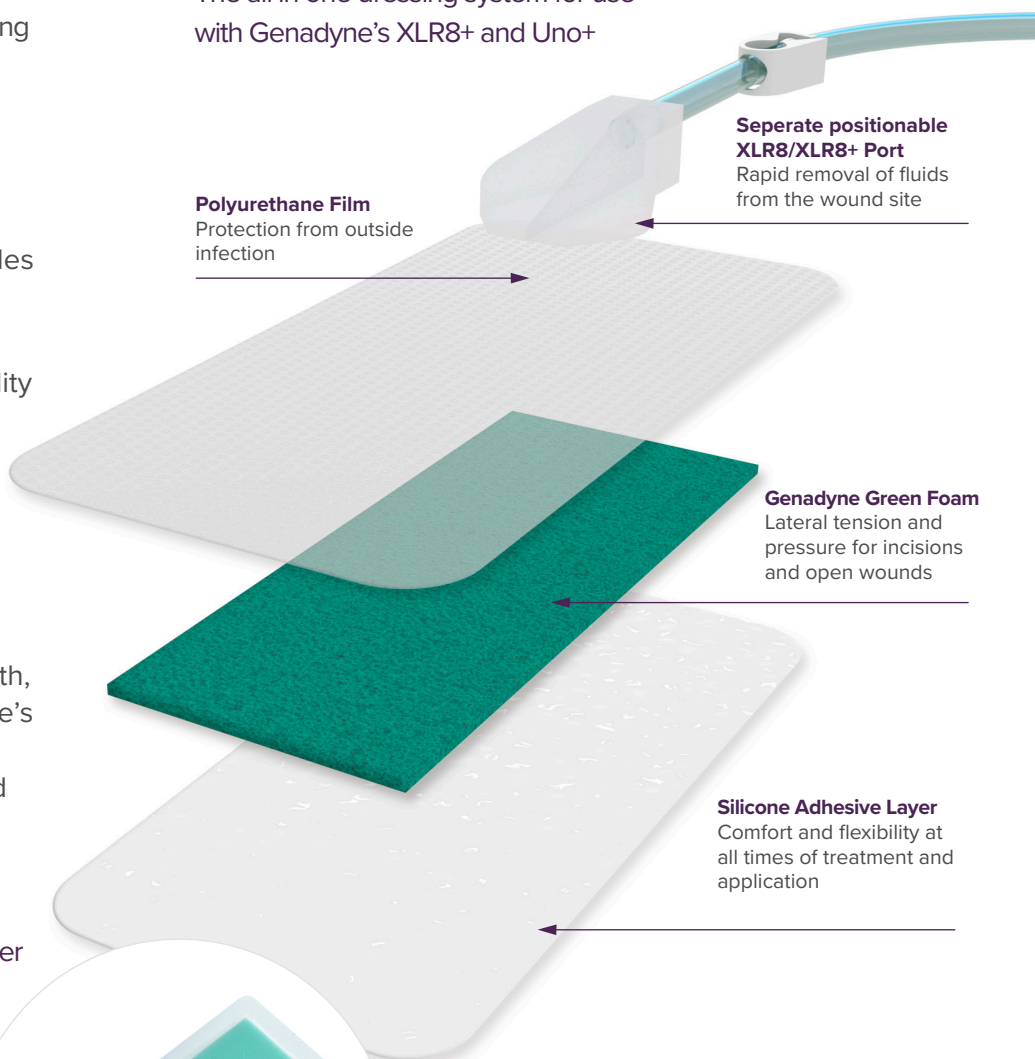

Separate positionable XLR8/XLR8+ Port
Rapid removal of fluids from the wound site

Polyurethane Film
Protection from outside infection

Genadyne Green Foam
Lateral tension and pressure for incisions and open wounds

Silicone Adhesive Layer
Comfort and flexibility at all times of treatment and application

Phone: 732-367-5533 www.handrhealthcare.com

the power to heal

DMR-13-093 Rev A